



IRONEDGE
MASSAGE GUN
PRO²

Instruction Manual

Introduction to products



Product name: Massage Gun Pro 2

Brand: IRONEDGE

Weight: 0.98kg

Dimension: 245*78*255mm

Input voltage: 110-240V, 50/60Hz

Plug type: DC, 16.8/1A

Massage mode: Warm compress, constant speed, undulating and climbing mode

Speed: 5 speed settings

Power: ≤60W

Rotate Speed: 1600-3200RPM/Min

Main material: PC+ABS

Power supply: 16.8V/2600mAh Rechargeable


Lithium-ion battery



Product structure




Battery and massage head installation method

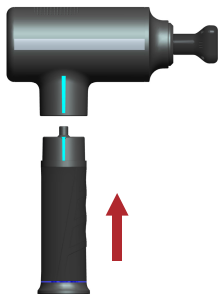


Battery installation method

- ★ Make the reading of main body align to the fluting of battery case

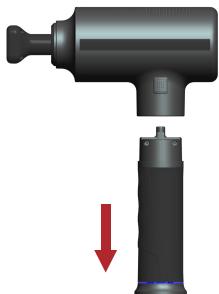

Battery Removal Method

- ★ Press the battery remove button, Pull the battery out

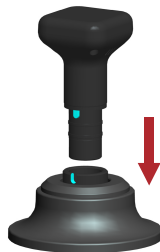

Heating head installation method

- ★ Make the reading of the heating head align to the fluting of the machine case



'U' shape head installation method


- ★ Make the reading of the 'U' shape head align to the fluting of the machine case



Remark:1.Other massage heads is not installed limit,direct vertical press into place.
2.Disassembly method of massage head:Pull them out vertically


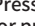
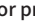
Product instructions





1、Turn on the main power switch at the bottom of the handle, Power indicator displays the current power indicator. (More than 40% of the power is displayed blue light; Less than 30% of the power is displayed yellow light; less than 10% of the power is displayed red light)

2、When the machine is in standby, press  to start massage, Then adjust the massage mode and speed according to your preference.

3、If you need to use 'Warm compress Mode'. Turn off the main power switch first, Pull out the existing massage head, Then take out the hot compress massage head and aim at the limit position. Turn on the main power switch, Press  to start Warm compress static massage mode, you can press  button to adjust the warm compress mode with the other three massage mode and speed according to your preference.

4、When the machine is working, Press  to pause the machine, When the machine is paused, Press  to start 'Warm compress static massage mode or press  to Return to previous working mode

5、The machine will turn off the mode display and speed display automatically when the machine is stop in 5 minutes. Press  button to start massage.

6、The machine is equipped with intelligent protection mode, The machine will automatically enter the sleep state When the motor works continuously for 20 minutes. It will take 1 minutes to restart. (Press  button to start massage)

Massage heads and Feature



Warm massage head

Temperature ranges from 42-48 °C, you can feel warm and comfortable, it benefits the blood circulation, it helps relax sore and stiff muscles and relieves the pain.



Round head

Used for hand, back, buttock, leg massage



Taper head

Used for deep tissue, such as joints, palm and plantar



Flat head

Used for relaxing and shaping each muscle part



“U” Shape head

Used for neck, spinal and achilles massage

Speed up the preheating effect before exercise



Activate the muscles

Keep your muscles active during exercise



Persistent state

Lactic acid is rapidly decomposed after exercise



Remove Fatigue

Shaking relaxes the myofascial membrane



Relieve Muscle Fatigue

Stretch the lines of the muscles



Fitness plasticity

Improves muscle balance



Auxiliary training

Deep release fatigue



Effective muscle relaxation

Relieve muscle tension



Relieve pain and stiffness

Rapid release of lactic acid



Accelerate muscle recovery

Instruction Manual



Main machine



Battery



Charger



Warm massage head



Taper head



Round head



Flat head



“U” Shape head

Maintain

If the following abnormal conditions occur during the use of this product, please refer to the solution:

Abnormal condition	Abnormal cause	Exception handling
Cannot start function	The battery may be low	Please charge
Automatic stop in use	1. The battery may be low	Please charge
	2.The motor has been working for 20 minutes	Restart in one minutes
The hot compress head temperature is not enough	Preheat the massage head for 1-2 minutes	Normal
The temperature of the hot compress massage head is too high	Hot compress massage mode is turned on without contact with the body to balance the temperature	Turn off the warm compress mode for 3 minutes before resuming normal use
The instrument is suspended for 5 minutes, and all the display lights on the main control board are off	Enter smart power saving mode	Press the start button to return to the previous working mode
Remark: If there are any other abnormal conditions that cannot be solved, please contact the customer service in time		

Matters needing attention

Disclaimer

- 1、Do not use this machine while charging
- 2、Do not insert or touch the gap between the massage head and the shell with your fingers or other objects during use, so as not to pinch.
- 3、Do not allow children to use alone.
- 4、This machine can not be self-disassembled or modified.
- 5、Due to the high power of this machine,it's recommended to use it for 15-20 minutes at a time to avoid physical discomfort.
- 6、If you find any abnormal phenomena or faults during use, please stop using them immediately and contact our company or dealer in time.
- 7、Be sure to use the original power adapter to charge the machine.
- 8、When the massage head is changed, it must be switched off and replaced.
- 9、Turn off the main power switch after each use of the instrument.
- 10、Do not rinse or soak with water or other liquids while cleaning. The best way is to use a wrung dry wet towel or medical alcohol cloth to wipe, to prevent water into the body.

Not suitable to use this device

Osteoporosis patients, people with heart problems and wearing electronic medical devices and pacemakers pregnant women or menstrual women, people with physical trauma or feel unwell by doctors, and are banned from this machine when the body is wet.